

# Indirect Calorimetry IMPC\_CAL\_003

## Purpose

Indirect calorimetry provides detailed information on the energy metabolism of mutant mice. Energy expenditure is evaluated through indirect calorimetry by measuring oxygen consumption with an open flow respirometric system. CO<sub>2</sub> and O<sub>2</sub> sensors measure the difference in CO<sub>2</sub> and O<sub>2</sub> concentrations in air volumes flowing through control or animal cages. The amount of oxygen consumed over a given period of time can thus be calculated, as far as the air flow through the cage is known. Data are expressed as ml O<sub>2</sub> h<sup>-1</sup>animal<sup>-1</sup>. The system also monitors CO<sub>2</sub> production, therefore, the respiratory exchange ratio (RER) and heat production can be calculated. An activity and food and water intake monitoring system can also be integrated into the set up in order to investigate circadian pattern and behaviour.

Ontological description: MP:0005266 - abnormal metabolism.

## Experimental Design

**Minimum number of animals :** 7M or 7F

**Age at test:** Week 11

Sexual dimorphism: In general, female mice have higher metabolism compared to males therefore statement is not entirely correct. However, genotype x sex interaction are rare therefore testing only males is acceptable.

It is essential that all phenotyping experimentation is conducted at the same time of day because physiological and biochemical parameters e.g. metabolic rate, body temperature and activity are subject to temporal rhythms. In the indirect calorimetry module standard measurements begin five hours before lights-off (lights off = T0) and are finished at T16 i.e. four hours after lights-on the next morning. Optional: Mice can be given one day of acclimation before the trial, and the trial can be continued for more than 21 hours.

## Equipment

1. Calorimetric system equipped with respirometer, feeder and water bottles
2. Ambulatory activity monitor (dependent on system specifications)
3. Food and water intake monitor
4. Computer with apparatus software installed

## Procedure

1. Optionally mice are allowed to acclimatise to the phenotyping room, to the calorimetry cage, food hoppers and drinking bottles 24 hours before testing.
2. Prepare and calibrate the calorimetric apparatus to confirm the accuracy of the gas sensors and flow meters. Specifically prior to each experiment:
  - a. Apply known volumes of CO<sub>2</sub> and O<sub>2</sub> to determine the sensitivity of the gas sensors and flow meters.
  - b. Run a complete calibration protocol according to the manufacturer's recommendations.
3. Provide each calorimetry cage with sufficient food and water for a period of ~24 (or 48) hours.
4. Weigh the mouse.
5. Place the mouse into a calorimetry cage with food and water available *ad libitum*.
6. Label the chamber with the corresponding subject identification and close it ensuring there is adequate air flow.
7. Initiate the calorimetric system for measurement:
  - a. Set up a new experiment in accordance with the manual (or load a file from a previous experimental setting).
  - b. Start recording measurements five hours before lights off for a total duration of 21 hours at minimum. Optional: 24 hours acclimation can be applied and the recording may continue for 48 hours.
  - c. The latency of CO<sub>2</sub> and O<sub>2</sub> activity transmitted and recorded is dependent on the number of chambers in use but will be logged periodically.
8. Generating a data report:
  - a. Upload all data from the experimentation including:
    - Gas analysis VO<sub>2</sub> and VCO<sub>2</sub> (ml/h/animal)
    - Heat production (kJ/h/animal)
    - Periodicity of measurements taken throughout experimentation (Figure 1)
    - Animal and the corresponding chamber that was used
    - The respiratory exchange ratio (RER) can be calculated using the VCO<sub>2</sub>/VO<sub>2</sub> ratio.
9. Activity parameters recorded will depend on the specification of calorimetric system used:
  - a. Ambulatory activity can be derived from the number beam splits during the session
  - b. Total activity can be derived from the number of fine movement (e.g. grooming behaviours) as well as ambulatory activity
  - c. An average of each of these parameters of activity is calculated hourly across the measurement period (between T-5 and T16).
  - d. Water and food intake (cumulative, hourly or total food and water intake, between T-5 and T16, will be computable depending on the calorimetric system used).

10. Remove each mouse from its chamber in turn at the end of the experimental session and record its weight. Return to their home cage.
11. Monitor the animals carefully to observe any abnormal behaviour(s). Ensure that food and water are available *ad libitum*.
12. Wash and wipe clean the chambers with warm water and dilute alcohol or appropriate disinfectant respectively.

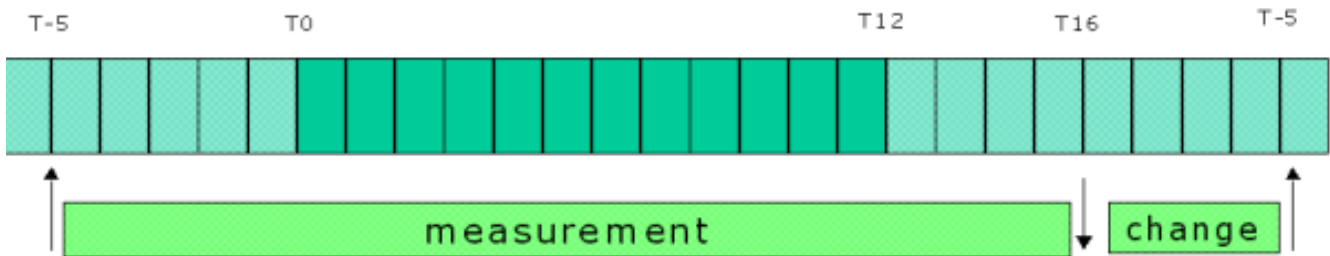


Figure 1. Daily workflow of calorimetric experimentation (Note: T0 designates start of dark cycle).

## Notes

The system requires periodic calibration of the gas sensors and flow meters to ensure precise measurements. The calibration procedure consists of the application of a gas of known composition and adjusting control knobs in the front of the oxygen and carbon dioxide sensors to obtain readings that reflect the contents of the calibration gas. It is recommended that the system be calibrated prior to the start of each experiment. The analyzers should not be shut down if not urgently required for maintenance. If this has to be done a warm up time of at least 90 minutes is required for the gas sensors for calibration (refer to manufacturer's manual). Calibrations and shut downs should be recorded in the laboratory journal.

Calorimetry test is to be performed before the ECG/ECHO test to avoid effects of hair removal on the calorimetry results.

The information about the date of the experiment, that is the date when the measurement is performed, is an important parameter which is to be submitted in the Experiment xml file (dateOfExperiment="2013-02-28").

## Data QC

1. Respiratory Exchange Rate (RER) is between 0.7-1.00
2. Mice show normal feeding and drinking behaviour
3. Mice show stable weight before and after calorimetry
4. Correct calibration of gases according to manufacturer's manual

## MetaData and examples

--	--

Metadata	Example
Time of dark cycle start	The starting time of the dark cycle. E.g. 7 pm
Time of dark cycle end	The ending time of the dark cycle. E.g. 7 am
Room temperature	The range of min. and max. temperatures of the calorimetry chamber. If the calorimetry machine doesn't register the chamber temperature, the room temperature in which the machine resides can be reported instead. E.g. 20.0-24.0 (C°). If the temperature is constant throughout the experiment, a single temperature can be submitted. E.g. 20.0 (C°). Do not submit the subtraction between the 2 values!
Acclimation to calorimetry cages	E.g. Yes/No.
Duration of test	Duration of the test without including the acclimation period. Can be a minimum of 21 hours or more, when acclimation is done. E.g. 21 (hours).
Equipment ID	ID of the machine used when more than 1 is used having same model and manufacturer. E.g. machine 1, machine 2, machine Minnie, machine Mickey Mouse, etc.
Equipment manufacturer	Manufacturer of the equipment. E.g. TSE Systems GmbH.
Equipment model	Model of the equipment. E.g. Labmaster CaloSys.
Experimenter ID	An ID of any format to be used coherently both inside the same procedure and for all procedures. E.g. Harw_001, or 1/2/3.
Date equipment last calibrated	Most recent date in which the equipment (or any part of) used in the procedure was subject to a calibration event.
Outer dimension of cage	(L, cm) x (W, cm) x (H, cm). E.g. 21.4 x 11.5 x 13.3

Inner dimension of cage	(L, cm) x (W, cm) x (H, cm). E. g. 21.4 x 11.5 x 13.3
Height from platform** to lid assembly	cm. E.g. 8.6
Available space for mouse	(L, cm) x (W, cm) x (H, cm). E.g. 20.3 x 10.4 x 8.6
Infrared beam setting on X axis	Orizontal activity: locomotion. N beams with L cm spacing, H cm above platform**. E.g. 16, 1.3, 3.0
Infrared beam setting on Y axis	Orizontal activity: locomotion. N beams with L cm spacing, H cm above platform**. E.g. 16, 1.3, 3.0
Infrared beam setting on Z* axis	Vertical activity: rearing/jumping. N beams with L cm spacing, H cm above platform**. E.g. 16, 1.3, 6.8
Beam strip placement on exterior of chamber****	Left/right, front/back, both
Lightbeam wavelength	nm, E.g. 950
Scanning rate	Hz, E.g. 100
Refractory period	s, E.g. 0.8
Presence of bedding into the cage	Yes/No
Igloo in cage	Yes/No

\*Z has same direction as X only is higher with respect to the platform

\*\* the floor where the mouse stand

\*\*\* This value is the "Height from platform\*\* to lid assembly"

## Parameters and Metadata

**Body weight before experiment** IMPC\_CAL\_001\_001 | v1.3

simpleParameter

**Req. Analysis:** false

**Req. Upload:** true

**Is Annotated:** false

**Unit Measured:** g

**Description:** body\_weight\_before\_experiment

---

## **Body weight after experiment** IMPC\_CAL\_002\_001 | v1.3

simpleParameter

**Req. Analysis:** false

**Req. Upload:** true

**Is Annotated:** false

**Unit Measured:** g

**Description:** body\_weight\_after\_experiment

---

## **Oxygen consumption** IMPC\_CAL\_003\_001 | v1.2

seriesParameter

**Req. Analysis:** false

**Req. Upload:** true

**Is Annotated:** true

**Unit Measured:** ml/h/animal

**Description:** oxygen\_consumption

**Increments:** Minimum 21

---

## Carbon dioxide production IMPC\_CAL\_004\_001 | v1.2

seriesParameter

Req. Analysis: false

Req. Upload: true

Is Annotated: true

Unit Measured: ml/h/animal

Description: carbon\_dioxide\_production

Increments: Minimum 21

---

## Heat production (metabolic rate) IMPC\_CAL\_005\_001 | v1.2

seriesParameter

Req. Analysis: false

Req. Upload: true

Is Annotated: true

Unit Measured: kJ/h/animal

Description: heat\_production\_metabolic\_rate\_

Increments: Minimum 21

---

## Ambulatory activity (no. of beam cuts) IMPC\_CAL\_006\_001 | v1.2

seriesParameter

Req. Analysis: false

Req. Upload: false

Is Annotated: true

Unit Measured: count/hour

**Description:** ambulatory\_activity\_no\_of\_beam\_cuts\_

**Increments:** Minimum 21

---

## **Total activity (no. of fine movement + no. of beam cuts) IM**

PC\_CAL\_007\_001 | v1.2

seriesParameter

**Req. Analysis:** false

**Req. Upload:** false

**Is Annotated:** true

**Unit Measured:** count/hour

**Description:** total\_activity\_no\_of\_fine\_movement\_no\_of\_beam\_cuts\_

**Increments:** Minimum 21

---

## **Total food intake IMPC\_CAL\_008\_001 | v1.2**

simpleParameter

**Req. Analysis:** false

**Req. Upload:** false

**Is Annotated:** true

**Unit Measured:** g

**Description:** total\_food\_intake

---

## **Cumulative food intake IMPC\_CAL\_009\_001 | v1.2**



seriesParameter

**Req. Analysis:** false

**Req. Upload:** false

**Is Annotated:** true

**Unit Measured:** g

**Description:** cumulative\_food\_intake

**Increments:** Minimum 21

---

## Acclimation to respirometry cages IMPC\_CAL\_012\_001 | v1.0

procedureMetadata

**Req. Analysis:** false

**Req. Upload:** true

**Is Annotated:** false

**Description:** acclimation\_to\_respirometry\_cages

**Options:** Yes, No,

---

## Duration of test IMPC\_CAL\_013\_001 | v1.0

procedureMetadata

**Req. Analysis:** false

**Req. Upload:** true

**Is Annotated:** false

**Unit Measured:** Hours

**Description:** duration\_of\_test

**Options:** 21, 48, 22, 23,

---

---

## Equipment ID IMPC\_CAL\_014\_001 | v1.0

procedureMetadata

**Req. Analysis:** false

**Req. Upload:** true

**Is Annotated:** false

**Description:** equipment\_name

---

## Equipment manufacturer IMPC\_CAL\_015\_001 | v1.1

procedureMetadata

**Req. Analysis:** true

**Req. Upload:** true

**Is Annotated:** false

**Description:** equipment\_manufacturer

**Options:** Columbus Instruments, Sable Systems, TSE Systems GmbH, O'hara Co. Ltd.,

---

## Equipment model IMPC\_CAL\_016\_001 | v1.2

procedureMetadata

**Req. Analysis:** true

**Req. Upload:** true

**Is Annotated:** false

**Description:** equipment

**Options:** Oxymax FAST, Oxymax/CLAMS, SAMPLEMAX, SM-MARS 8 channel Metabolic system, PhenoMaster/LabMaster CaloSys, (Dri)+Fed, Act.X,Y 16mice, PhenoMaster/Labmaster CaloSys, (Dri)+Fed, Act.X,Y 12Mice, FWI-3002 & IA-16M,

PhenoMaster/Labmaster CaloSys, (Dri)+Fed, Act.X,Y 20Mice,  
PhenoMaster/Labmaster CaloSys, (Dri)+Fed, Act.X,Y 32Mice,  
PhenoMaster/Labmaster CaloSys, (Dri)+Fed, Act.X,Y 12Mice - Gas system 994620 series,  
PhenoMaster/Labmaster CaloSys, (Dri)+Fed, Act.X,Y 8Mice,

---

## Respiratory Exchange Ratio IMPC\_CAL\_017\_001 | v1.2

simpleParameter

Req. Analysis: false

Req. Upload: false

Is Annotated: true

Description: respiratory\_exchange\_ratio

Derivation:

div(meanOfIncrements('IMPC\_CAL\_004\_001',21), meanOfIncrements('IMPC\_CAL\_003\_001',21))

---

## Experimenter ID IMPC\_CAL\_018\_001 | v1.0

procedureMetadata

Req. Analysis: false

Req. Upload: true

Is Annotated: false

---

## Date equipment last calibrated IMPC\_CAL\_019\_001 | v1.2

procedureMetadata

Req. Analysis: false

Req. Upload: false

Is Annotated: false

---

## Total water intake IMPC\_CAL\_021\_001 | v1.1

simpleParameter

Req. Analysis: false

Req. Upload: false

Is Annotated: true

Unit Measured: ml

---

## Cumulative water intake IMPC\_CAL\_022\_001 | v1.1

seriesParameter

Req. Analysis: false

Req. Upload: false

Is Annotated: true

Unit Measured: ml

Increments: Minimum 21

---

## Room temperature max IMPC\_CAL\_023\_001 | v1.4

procedureMetadata

Req. Analysis: false

Req. Upload: true

Is Annotated: false

Unit Measured: C

---

## Room temperature min IMPC\_CAL\_024\_001 | v1.2

procedureMetadata

Req. Analysis: false

Req. Upload: true

Is Annotated: false

Unit Measured: C

---

## Chamber temperature min IMPC\_CAL\_025\_001 | v1.1

procedureMetadata

Req. Analysis: false

Req. Upload: false

Is Annotated: false

Unit Measured: C

---

## Chamber temperature max IMPC\_CAL\_026\_001 | v1.1

procedureMetadata

Req. Analysis: false

Req. Upload: false

Is Annotated: false

Unit Measured: C

---

## Outer dimension of cage IMPC\_CAL\_027\_001 | v1.0

procedureMetadata

Req. Analysis: true

Req. Upload: false

Is Annotated: false

**Unit Measured:** cm

**Description:** (L, cm) x (W, cm) x (H, cm)

---

## **Inner dimension of cage** IMPC\_CAL\_028\_001 | v1.0

procedureMetadata

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false

**Unit Measured:** cm

**Description:** (L, cm) x (W, cm) x (H, cm)

---

## **Height from platform to lid assembly** IMPC\_CAL\_029\_001 | v1.0

procedureMetadata

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false

**Unit Measured:** cm

---

## **Available space for mouse** IMPC\_CAL\_030\_001 | v1.0

procedureMetadata

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false

**Unit Measured:** cm

**Description:** (L, cm) x (W, cm) x (H, cm)

---

## **Infrared beam setting on X axis** IMPC\_CAL\_031\_001 | v1.0

[procedureMetadata](#)

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false

**Unit Measured:** cm

**Description:** Horizontal activity: locomotion. N beams with L cm spacing, H cm above platform

---

## **Infrared beam setting on Y axis** IMPC\_CAL\_032\_001 | v1.0

[procedureMetadata](#)

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false

**Unit Measured:** cm

**Description:** Horizontal activity: locomotion. N beams with L cm spacing, H cm above platform

---

## **Infrared beam setting on Z axis** IMPC\_CAL\_033\_001 | v1.0

[procedureMetadata](#)

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false

**Unit Measured:** cm

**Description:**

Vertical activity: rearing/jumping. N beams with L cm spacing, H cm above platform.

---

**Beam strip placement on exterior of chamber** IMPC\_CAL\_034\_0

01 | v1.0

[procedureMetadata](#)

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false

**Options:** Left/right, Front/back, Both,

---

**Lightbeam wavelength** IMPC\_CAL\_035\_001 | v1.0

[procedureMetadata](#)

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false

**Unit Measured:** nm

---

**Scanning rate** IMPC\_CAL\_036\_001 | v1.1

[procedureMetadata](#)

**Req. Analysis:** true

**Req. Upload:** false

**Is Annotated:** false



Unit Measured: Hz

---

## Refractory period IMPC\_CAL\_037\_001 | v1.0

procedureMetadata

Req. Analysis: true

Req. Upload: false

Is Annotated: false

Unit Measured: s

---

## Presence of bedding into the cage IMPC\_CAL\_038\_001 | v1.0

procedureMetadata

Req. Analysis: true

Req. Upload: false

Is Annotated: false

Options: Yes, No,

---

## Igloo in cage IMPC\_CAL\_039\_001 | v1.0

procedureMetadata

Req. Analysis: true

Req. Upload: false

Is Annotated: false

Options: Yes, No,

---

## Calibration method IMPC\_CAL\_040\_001 | v1.0

procedureMetadata

Req. Analysis: true

Req. Upload: false

Is Annotated: false

Options: Reference range method, Gas bottle method,

---

## Acclimation period IMPC\_CAL\_041\_001 | v1.0

procedureMetadata

Req. Analysis: false

Req. Upload: false

Is Annotated: false

Unit Measured: Hours

---

## Time of dark cycle start IMPC\_CAL\_010\_003 | v3.0

procedureMetadata

Req. Analysis: false

Req. Upload: true

Is Annotated: false

### Description:

The starting datetime of the dark cycle. E.g. 2013-02-28T06:25:00-08:00+01 or 2013-02-28T06:25:00Z

---

## Time of dark cycle end IMPC\_CAL\_020\_003 | v3.1

procedureMetadata

**Req. Analysis:** false

**Req. Upload:** true

**Is Annotated:** false

**Description:**

The ending time of the dark cycle. E.g. 2013-02-28T06:25:00-08:00+01 or 2013-02-28T06:25:00Z

---

**Respiratory Exchange Ratio series** IMPC\_CAL\_042\_001 | v1.0

seriesParameter

**Req. Analysis:** false

**Req. Upload:** false

**Is Annotated:** false

**Increments:** Minimum 21

---